



ALL-IN-ONE – HSW MED-CAM SYSTEM

- Ultra-slim HSW Video Nasolaryngoscope
- High-quality HD camera head
- All-purpose system for ENT applications

All-in-one – All you need

This all-in-one HSW Med-cam system has been developed in close cooperation with clinical users in order to provide maximum benefits for both patients and users. The ultra-slim HSW Video Nasalaryngoscope provides an

extremely gentle insertion, meaning minimal discomfort and time in procedure. The high-quality HD camera head offers all possibilities and enables optimal image quality.


CMOS Video Nasalaryngoscope – Optimal cost/performance ratio



Integrated LED
for optimal operating expense



Ergonomic design
for effortless use

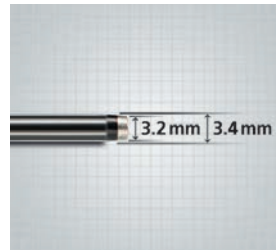

Ultra-slim-line
for gentle insertion



Precise motion
for perfect control

CMOS HD camera head – All possibilities



HSW Video Nasalaryngoscope with integrated LED illumination

Article no.	830 001 1825
Working length	approx. 300 mm
Diameter distal end	3.2 mm
Diameter insertion tube	3.4 mm
Direction of view	0°
Field of view	90°
Bending range of distal tip	130° up/down

HSW Video Nasalaryngoscope, light guide version

Article no.	830 001 2731
-------------	--------------

HSW camera head

Article no.	830 001 3118
-------------	--------------

HSW coupler

Article no.	830 001 8854
-------------	--------------

HSW Med-cam control unit

Article no.	830 001 2748
Input voltage range	100–240 VAC, 50/60 Hz
Output	1 x USB for image storage, DVI